SLIT LAMP SL-500/SL-203/SL-102

SL-500 SPECIFICAT	TON		
MICROSCOPE	Туре	Binoculer, Galilean-type, 5step drums rotation	
	Conversion	12 degrees	
	Eyepieces	12.5x High-Eyepoint	
	Filter	Yellow filter	
	Total Magnification	6x, 10x, 16x, 25x, 40x	
	Field of View	φ37, φ23, φ14.6, φ9.5, φ5.9mm	
	Pupillary Distance	52mm to 80mm	
	Diopter Adjustment Range	-7D to +7D	
Base	Side Shift	100mm	
	Depth Shift	100mm	
	Height Shift	29mm	
Chin-rest	Height Shift	80mm	
	Fixation Target	Red LED	
	Slit Width	0 to 15mm, continuous	
ILLUMINATION SYSTEM	Slit Height	0 to 15mm, continuous	
	Diaphragm Sizes	φ0.2, φ1, φ3, φ5, φ8, φ15mm	
	Rotation	0 to 180 degrees	
	Filters	Cobalt blue, 50%ND, Heat absorption, Green(Red-Free)	
	Light Source (Bulb)	12V 30W Halogen lamp	
Power Consumption		60VA	
Size	Weight	23kg	
	Dimensions	550mm(W) × 430mm(D) × 780mm(H)	

SL-500	
Standard	Accessories

- Operation Manual
- Power Cable
- Dust Cover
- Chin-rest Paper
- Halogen Lamp
- Mirror
- Test Bar
- Light Guide
- Hand Rest

SL-500 Optional Parts

- Beam Splitter
- C-mount
- CCD Camera
- Observation Tube
- Hruby lens attachment ■ Applanetion Tonometer

SL-203/102 SPECIFICATION		SL-203	SL-102	
MICROSCOPE	Туре	Galilean-type	Greenough-type	
	Magnification Changing Method	5step drums rotation	2-position lever	
	Eyepieces	15x	10x, 16x	
	Total Magnification	6x, 10x, 16x, 25x, 40x	10x, 16x, 25.6x	
	Field of View	6~37mm	10~18mm	
	Pupillary Distance	50~80mm	55~80mm	
	Diopter Adjustment Range	-7D to +7D	-8D to +8D	
Base	Side Shift	100mm		
	Depth Shift	100mm		
	Height Shift	29mm		
Chin-rest	Height Shift	80mm		
	Fixation Target	Red LED		
ILLUMINATION SYSTEM	Slit Width	0 to 15mm, continuous		
	Slit Height	0 to 15mm, continuous		
	Diaphragm Sizes	φ0.2, φ1, φ3, φ5, φ8, φ15mm		
	Rotation	0 to 180 degrees		
	Filters	Cobalt blue, 50%ND, Heat absorption, Green		
	Light Source (Bulb)	12V 30W Halogen lamp		
Power Consumption		40VA		
Size	Weight	18.5 kg	18 kg	
	Dimensions	550mm(W) x 400mm(D) x 782mm(H)		

SL-203 / SL-102 **Standard Accessories**

- Operation Manual
- Power Cable
- Dust Cover
- Chin-rest Paper ■ Halogen Lamp
- Mirror
- Test Bar

SL-203 Optional Parts

- Beam Splitter
- C-mount
- CCD Camera ■ Observation Tube
- Eye piece W/Yellow Filter
- Hruby lens attachment
- Applanetion Tonometer

SL-102 Optional Parts

- Hruby lens attachment
- Applanetion Tonometer



Ohira Co.,Itd.

629-961 Niiborishinden, Minamiuonuma, Niigata, 949-7135, Japan

Distributed by



MEC Sales Division 2-8-4, Kandatsukasa-machi, Chiyoda-ku Tokyo, 101-0048, Japan

TEL 81-3-3256-7701 FAX 81-3-3256-7702

Design and specifications are subject to change without notice.

E-mail: eye@rexxam.co.jp URL: http://www.rexxam.co.jp URL: http://www.shin-nippon.jp




SHIN-NIPPON by **Rexxam**

SL-500

Upgraded illumination system! Brighter and more user-friendly.

SL-500

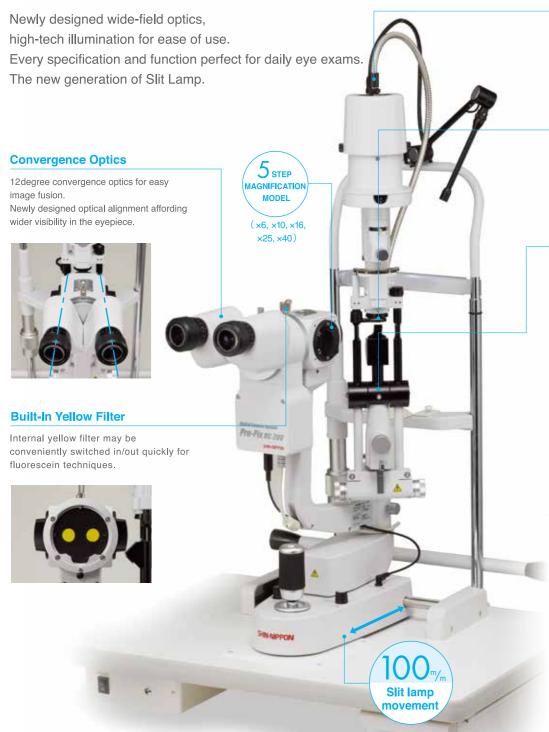

Illumination Source

The new improved illumination pathway is created to increase overall luminance.
The luminance output of the unit's 30W Halogen Lamp is comparable to our existing models with a 50W Halogen Lamp.

Expanded Slit Width & Height

A wider slit width & longer slit height capability are attractive features when combined with our new wide field optics. Slit width:0-15mm
Slit height:0-15 mm

Background illumination/Diffuser

Flexible background illumination and a diffuser that can be easily deployed, e.g. when shooting.



Improved Working Distance

Forward and back slit lamp movement is increase up to 100mm allowing retinal exams with diagnostic and therapeutic fundus lenses.



100mm

Wide Field & High Eye point Eyepiece Design



With our newly created wide field optics, we have successfully increased the actual field of view. Wider fields of view as 37, 23, 14.6, 9.5 & 5.9 mm are achieved. With the new high-point eyepiece feature, eye exams are performed in a more natural, comfortable manner.

Centralized Control



Joystick with camera shutter control button, brightness control knob and locking control knob, routinely used functions all conveniently Localized in the joystick area for efficient, confident slit lamp control.

SL-203


The acclaimed SL-203 is now upgraded (a stereotype eye lens tube). Equipped with the same illumination system as the SL-500.

The main unit can now slide up to 100 mm,

which makes the device brighter and more user-friendly.

Customizable in various ways to meet your needs,

e.g. allowing a digital camera to be added to the basic function.

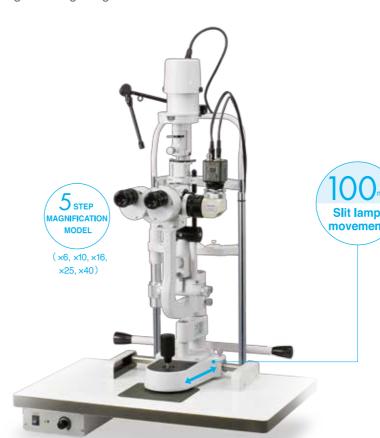





SL-102

The simple and easily operable SL-102
Also with upgraded illumination system and expanded slide travel of the main unit.
A basic model for all users,

allowing you to observe the eye quickly and correctly, even when multiple examinations are required.


DIGITAL CAMERA SYSTEM

Digital Camera System "Pro-Pix DC-200" is available as an option for Model: SL-500 and SL-203.

Database Software - Standard with Pro-Pix

Software which can be installed into your PC for Photo and Video of patient examinations. Stored images may be edited.



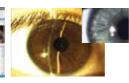

Power Cable Not Required

Pro-Pix Camera is powered by the USB port of your PC, helps to reduce additional cable clutter.



 Camera
 1/2inch CMOS camera 1600 x 1200 pixel

 Power source
 USB 2.0 5V 4A

